MEDICAL & BEAUTY



OCEAN

MB/PC280
INSTRUCTION MANUAL

WARNINGS

Safety during installation

• Before proceeding with the installation and use of the chair, carefully read the instruction manual, making sure to follow all instructions and warnings provided therein. It contains important safety information regarding the installation, use, and maintenance of the product!

• After removing the packaging, check that the contents are intact and complete. If you have any doubts, do not use the product and contact our customer service.

• The Ocean chair must be placed on a flat, level surface capable of supporting the weight of the chair when loaded, taking into account the maximum allowable load.

• Warning! For a correct and safe installation, it is recommended to move the Ocean chair as shown in Fig. 2

Safety during use

• The appliance must be used exclusively for the purpose for which it was expressly designed. Any other use is to be considered improper and therefore hazardous. The manufacturer cannot be held responsible for any damage resulting from improper, erroneous, or unreasonable use, which will void all warranty rights.

• Use is restricted exclusively to personnel performing the treatment (podiatrist, beautician).

• Children must not be left unattended near the appliance.

Safety during maintenance

• Cleaning and maintenance must be performed only by personnel qualified for such activities and must not be carried out by persons unfamiliar with its operation.

Transport / Storage / Operating requirements

Transport must take place exclusively inside the original packaging.

Transport and storage temperature: from -20 to +50C

Operating temperature: from +10 to +40C

Relative air humidity: max. 90%, non-condensing

Avoid sudden temperature variations to prevent condensation.

INTENDED USE

The Ocean chair has been designed to offer maximum comfort during podiatry and pedicure treatments. By rotating 90°, it ensures ease of use for both the client and the operator. The chair rests on a rectangular stainless steel structure housing a white ceramic pedicure basin, complete with a mixer and floor drainage system. The Ocean is also equipped with an upholstered footrest that is adjustable in height and fore-and-aft.

COMPONENTS

Fig. 1 shows the main components of the Ocean podiatry treatment chair.

Legend:

1. Stainless steel base
2. White ceramic pedicure basin
3. Hardware: mixer tap with swivel spout
4. Adjustable footrest
5. Swivel chair (90°)


Fig. 1







Fig. 2

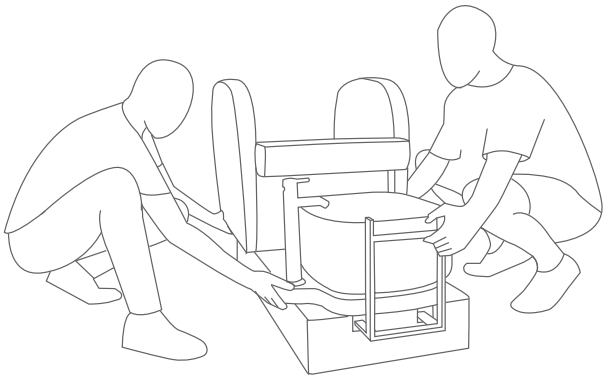


To proceed with the installation:

1. Remove the rear panel of the chair
2. Remove the rear and front stainless steel panels of the base
3. Position the semi-rotation stop on the desired side
4. Adjust the adjustable feet, if necessary
5. Mark the floor, drill, and anchor the structure using the four adjustable corner brackets
6. Connect the drainage and water supply, see Fig. 4
7. Check the watertightness of the water connection and drainage
8. Reinstall the panels

Fig. 4

A = COLD WATER SUPPLY 3/8”
B = DRAINAGE Ø40
C = HOT WATER SUPPLY 3/8”

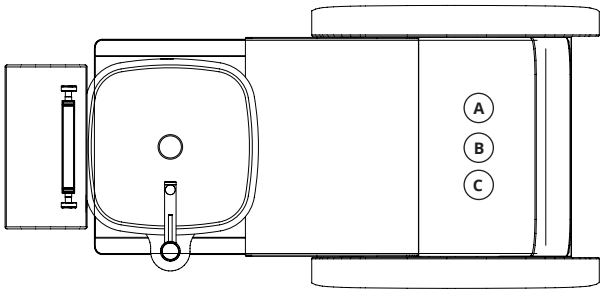

USE

To proceed with use:

1. Rotate the chair to allow the client to sit down
2. Adjust the height and depth of the footrest cushion using the appropriate handwheels
3. Proceed with the treatment
4. The drain waste plug is equipped with a click-clack system
5. Sanitize the basin after each use

MAINTENANCE

Cleaning and maintenance

Clean the chair upholstery using only lukewarm water, mild soap, and a soft cloth.

Stains must be removed immediately to prevent them from penetrating the fabric.

Clean the chair frame with a slightly damp cloth.

Do not rub disinfectants onto the upholstery.

Do not use chemicals or aggressive detergents that could damage the upholstery.

Do not use substances containing oil or grease to clean the upholstery.

The upholstery is not resistant to chemicals or dry cleaning.

Repair and maintenance

Periodically check that the water fittings are working properly.

The chair does not require any special maintenance if used correctly and in compliance with the manufacturer’s instructions and guidelines. However, in the event of any problems, please contact our customer service.

All technical interventions must be performed by qualified professionals.

Disposal

The various packaging components, the chair, the upholstery, and their relative accessories must be disposed of in compliance with the regulations established by local authorities.

EQUIPMENT FEATURES

MEDICAL & BEAUTY
100% MADE IN ITALY

Costruttore	Salon Habitat (Italia) s.r.l. via Enzo Lazzaretti 12 42122 Reggio Emilia www.medicalandbeauty.com info@medicalandbeauty.com
Modello	Ocean MB/PC280
Funzione	Support equipment for operators in aesthetic and podiatric treatments.
Anno di costruzione	2026
Massimo carico ammesso alla poltrona	130 Kg

MEDICAL & BEAUTY

100% MADE IN ITALY

Salon Habitat (Italia) s.r.l. via Lazzaretti 12 42122 Reggio Emilia
Sales Italia 0522 276740 International Sales 0522 276755 Fax 0522 276767
info@medicalandbeauty.com • www.medicalandbeauty.com